
**Jet dentaire
Oral Irrigator
Monddouche
Irrigador dental**

971037

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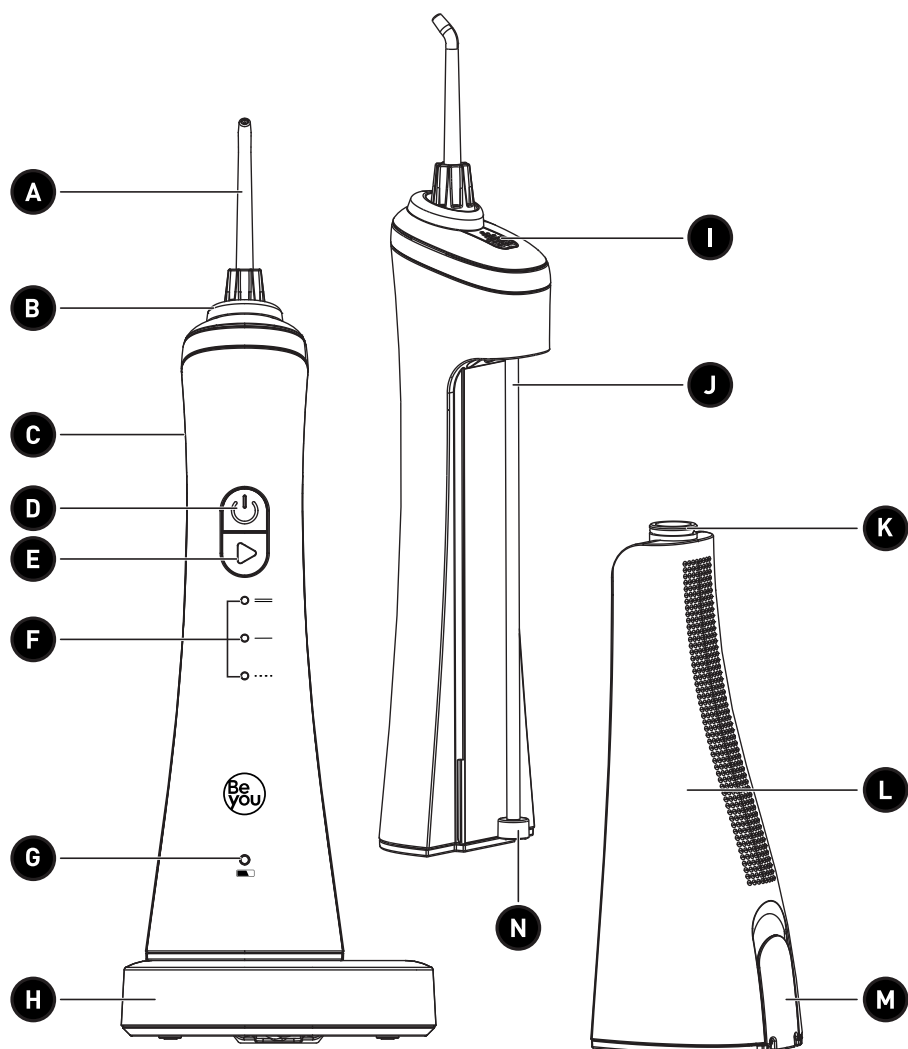Thank you!

Thank you for choosing this BEYOU product.
BEYOU products are chosen, tested, and recommended
by ELECTRO DEPOT, so you can be sure you are
getting a top-quality, easy-to-use product that won't
let you down.

We're confident your new device will be a pleasure
to use every day!

Welcome to ELECTRO DEPOT.

ELECTRO DEPOT

Table of Contents

A

Product overview

Product description

B

Product usage

Operating instructions

C

Cleaning and maintenance

Cleaning and maintenance
Troubleshooting
Disposing of the unit
Removing the rechargeable battery

Product description

A Jet tip

B Jet tip removal device

C Unit body

D Power button

E Mode button

F Mode indicator

G Charge indicator

H Wireless charger

I Jet tip removal button

J Suction tube

K Waterproof seal


L Water tank

M Water tank flap


N Filter

Operating instructions

1. Make sure your hands are dry before use, then plug the unit into a power source.

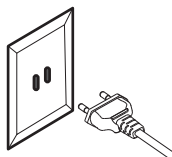
2. Insert the unit's body into the charger and the charging indicator will flash. It will continue to flash during the charging process. When the oral irrigator is fully charged, the charging indicator will become steadily lit.
It takes around 24 hours to fully charge the oral irrigator.

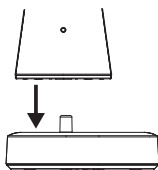
**NOTE:**

Should the charging indicator flash before or after use, this means that the oral irrigator needs to be charged.

3. Remove the oral irrigator from the charger, then insert the jet tip into the body of the unit.
To remove the jet tip, press the jet tip removal button.

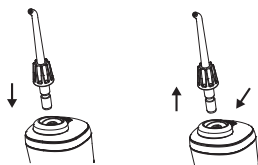4. The water tank can be filled in two different ways:
 - 4.1 Open the water tank valve and hold the water tank horizontally to fill it with water. Once filled, close the water tank valve.

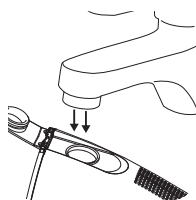- 4.2 Remove the water tank and hold it vertically to fill it with water. Then, reposition the water tank.


5. Press the "Mode" button to select the desired mode. There are three modes: "Normal", "Soft" and "Pulse" (see below). The "Mode" button can be pressed to switch between modes during use.
- A. **==** Normal: Normal teeth cleaning function.
 - B. **—** Soft: For sensitive teeth, first-time oral irrigator users, or if only a gentle clean is required.
 - C. **....** Pulse: To massage the gums, promote blood circulation and strengthen teeth.

**NOTE:**

When using this oral irrigator for the first time, please use the soft mode for the first week. Once you have become accustomed to the unit, you can use whichever mode you like.


6. Hold the oral irrigator vertically and insert the jet tip into your mouth; aim it at the gums or teeth and open your mouth slightly to let the water flow out.

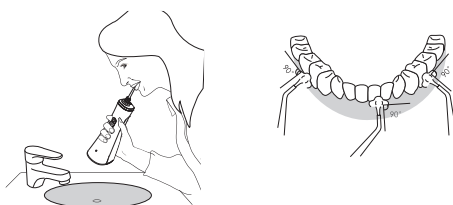7. Press the power button to switch on the oral irrigator. The following instructions should be taken into account during use.

A. The water should flow vertically over the gums or teeth.

B. Slowly move the oral irrigator over your teeth in a normal cleaning order. Do not spray water directly along the periodontal pocket.

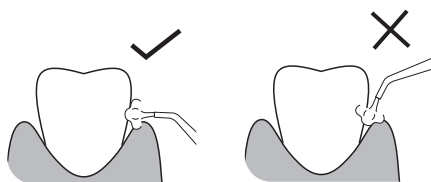
8. After use, press the power button to switch off the unit. The unit will automatically shut off after 2 minutes of use. Should you need to continue using the unit, press the power button again.

9. After use, please follow the instructions below to drain the residual water:

A. Open the water tank valve and empty the tank.


B. Press the power button; the oral irrigator will start to automatically drain the residual water.

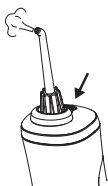C. Press the power button again to switch off the oral irrigator, then close the tank valve.

D. Wipe down the oral irrigator with a dry cloth.

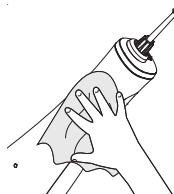**NOTE:**

To prevent bacteria proliferating, residual water should not be left inside the oral irrigator. Therefore, in the event of prolonged non-use, please thoroughly clean, dry and put away the oral irrigator.

Cleaning and maintenance

Please clean the oral irrigator with clean water or neutral detergent; never use corrosive or abrasive cleaners as these could damage the oral irrigator. Furthermore, do not wash it with water hotter than 50°C.

Cleaning the unit's body

- Please clean the unit's body with a dry or damp cloth.
- Do not immerse the unit's body in water when cleaning it.

Cleaning the water tank

- Wipe the water tank with a dry or damp cloth or detach it to rinse it with water.
- Fully drain the water tank in the event of prolonged non-use.

Cleaning the jet tip

- Wipe the jet tip with a dry or damp cloth or detach it to rinse it with water.
- Do not bend, pull or deform the jet tip.
- Dentists recommend replacing the jet tip every 6 months.

Cleaning the charger

Wipe the metal part of the plug with a dry cloth every 6 months to prevent the build-up of dust, which could, in turn, lead to contact issues.

Troubleshooting

Problem	Possible Cause	Solution
The oral irrigator is not working.	The oral irrigator is being used for the first time or has not been used for more than 3 months.	Please charge the oral irrigator for more than 24 hours.
The water pressure is low during use.	Something is blocking the jet tip. The jet tip is deformed.	Replace the jet tip.
The oral irrigator only works for a short time after being fully charged.	- The filter is clogged. - The battery is approaching the end of its life. - The oral irrigator is not charging properly.	- Clean the filter. - Please contact the after-sales service. - Please charge the oral irrigator for more than 24 hours before using it.
The charging indicator does not light up when charging the unit.	- The mains plug has not been connected properly. - Circuit malfunction.	- Please insert the mains plug correctly. - Please contact the after-sales service.


Disposing of the unit

- The oral irrigator is fitted with a recyclable rechargeable battery. When the oral irrigator comes to the end of its service life, remove the battery and recycle or dispose of it according to your local waste disposal regulations.
- Please consult your local council should you have any waste-related questions.


Removing the rechargeable battery

Only disassemble the oral irrigator before disposing of it, as this will otherwise damage the unit. Please unplug the unit before removing the battery.


- Make sure that the battery is dead before removing it. To fully discharge the battery, press the power button to start the oral irrigator and let it run until it stops automatically. Repeat the operation several times until the battery dies.
- Remove the water tank, then unscrew the two screws on the bottom cover with a Phillips screwdriver.

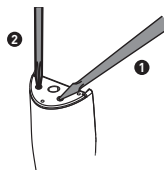3. Remove the two screws on the top cover with a Phillips screwdriver, then detach the top cover and remove the unit's central section.
4. Remove the screws located on the bracket of the unit's central section, by means of a Phillips screwdriver.

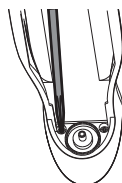5. Cut the battery connection wire with scissors, then remove the battery.

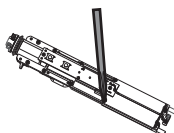6. Properly recycle or dispose of the battery. Please consult your local council for information regarding recycling the other oral irrigator components.

**NOTE:**

Do not short-circuit the battery and never turn the oral irrigator back on after removing the battery.

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